




## PneuGel™ Tennis Elbow Support

Single Patient Use Only

Rx Only

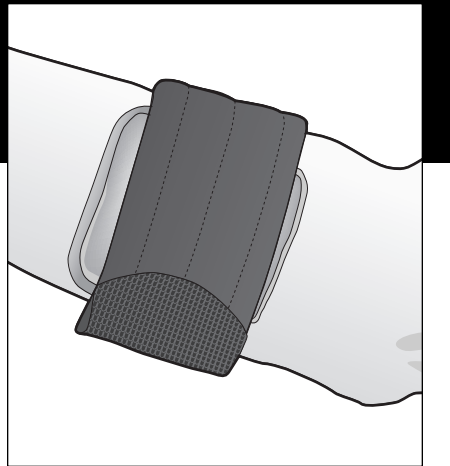

### IMPORTANT INFORMATION

Please read all instructions, warnings, and precautions before use.

### INTENDED USE



Intended for lateral and medial epicondylitis

### PRECAUTIONS



Consult your clinician immediately if you experience sensation changes, unusual reactions, swelling or increased pain while using this product. Use caution when fitting this product on patients with reduced sensitivity. Patients with circulation problems should not over tighten the strap. Product should be worn during activity only. Remove product when in a resting position or sleeping.

### DIRECTIONS FOR USE

1. With the gel bladder affixed to the strap, wrap the strap around the forearm about 1"-2" below the elbow.
2. Make sure the gel bladder is placed over the forearm muscle area and pull strap through d-ring closure. Secure with hook.

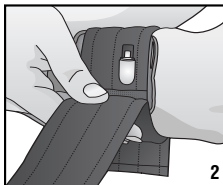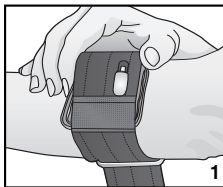

3. Inflate the gel bladder to desired pressure using the clear tip of the pump. If you over inflate simply insert the black tip of the pump to release the proper amount of air. **WARNING:** Over inflation could cause numbness and reduce circulation.

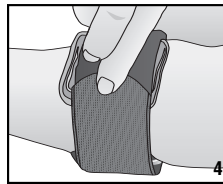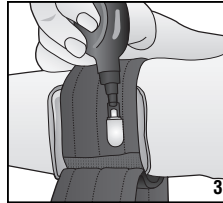

4. Use the remaining part of the strap to cover the valve and provide a secure fit. **DO NOT OVERTIGHTEN! SLIGHT PRESSURE IS SUFFICIENT.**

**COLD THERAPY INSTRUCTIONS:** The removable gel bladder may be placed in the freezer to provide cold therapy. The gel bladder can be placed in the freezer overnight, but it is generally ready for use within two hours in the freezer. Consult a physician prior to using cold therapy. **WARNING! THIS PRODUCT IS NOT DESIGNED FOR WARM THERAPY. DO NOT PLACE IN A MICROWAVE OVEN.**

### CLEANING INSTRUCTIONS


Remove gel bladder and hand wash strap with a mild detergent. Air Dry.



## PneuGel™ Tennis Elbow Support

Single Patient Use Only

Rx Only

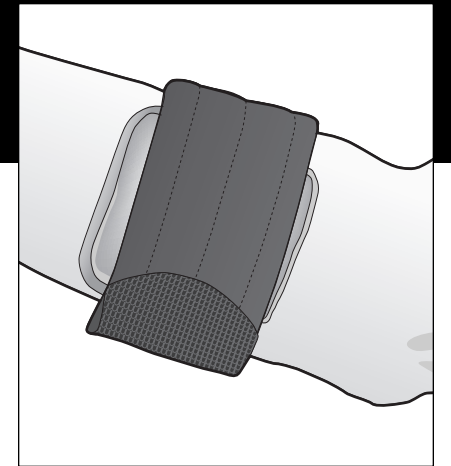

### IMPORTANT INFORMATION

Please read all instructions, warnings, and precautions before use.

### INTENDED USE



Intended for lateral and medial epicondylitis

### PRECAUTIONS



Consult your clinician immediately if you experience sensation changes, unusual reactions, swelling or increased pain while using this product. Use caution when fitting this product on patients with reduced sensitivity. Patients with circulation problems should not over tighten the strap. Product should be worn during activity only. Remove product when in a resting position or sleeping.

### DIRECTIONS FOR USE

1. With the gel bladder affixed to the strap, wrap the strap around the forearm about 1"-2" below the elbow.
2. Make sure the gel bladder is placed over the forearm muscle area and pull strap through d-ring closure. Secure with hook.

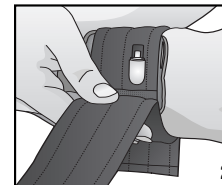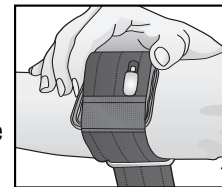

3. Inflate the gel bladder to desired pressure using the clear tip of the pump. If you over inflate simply insert the black tip of the pump to release the proper amount of air. **WARNING:** Over inflation could cause numbness and reduce circulation.

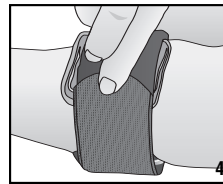
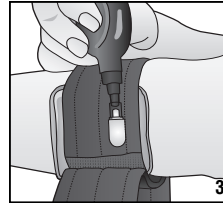


4. Use the remaining part of the strap to cover the valve and provide a secure fit. **DO NOT OVERTIGHTEN! SLIGHT PRESSURE IS SUFFICIENT.**

**COLD THERAPY INSTRUCTIONS:** The removable gel bladder may be placed in the freezer to provide cold therapy. The gel bladder can be placed in the freezer overnight, but it is generally ready for use within two hours in the freezer. Consult a physician prior to using cold therapy. **WARNING! THIS PRODUCT IS NOT DESIGNED FOR WARM THERAPY. DO NOT PLACE IN A MICROWAVE OVEN.**

### CLEANING INSTRUCTIONS

Remove gel bladder and hand wash strap with a mild detergent. Air Dry.

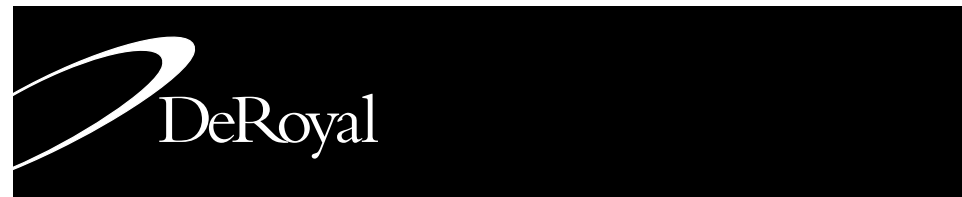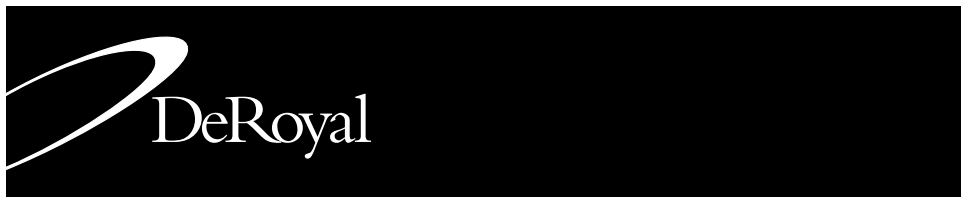

**Manufacturer:**  
Sroufe Healthcare Products Inc.  
P.O. Box 347  
Ligonier, IN 46767

**Distributor:**  
DeRoyal  
200 DeBusk Lane, Powell, TN 37849 U.S.A.  
888.938.7828  
865.938.7828

©2004 DeRoyal. All Rights Reserved.  
Part #0-1805  
Revised 9/04



**Manufacturer:**  
Sroufe Healthcare Products Inc.  
P.O. Box 347  
Ligonier, IN 46767

**Distributor:**  
DeRoyal  
200 DeBusk Lane, Powell, TN 37849 U.S.A.  
888.938.7828  
865.938.7828

©2004 DeRoyal. All Rights Reserved.  
Part #0-1805  
Revised 9/04