


Your Front-Line Defense for Optimum Fluid Control



Suction Mats & Anti-Fatigue Mats

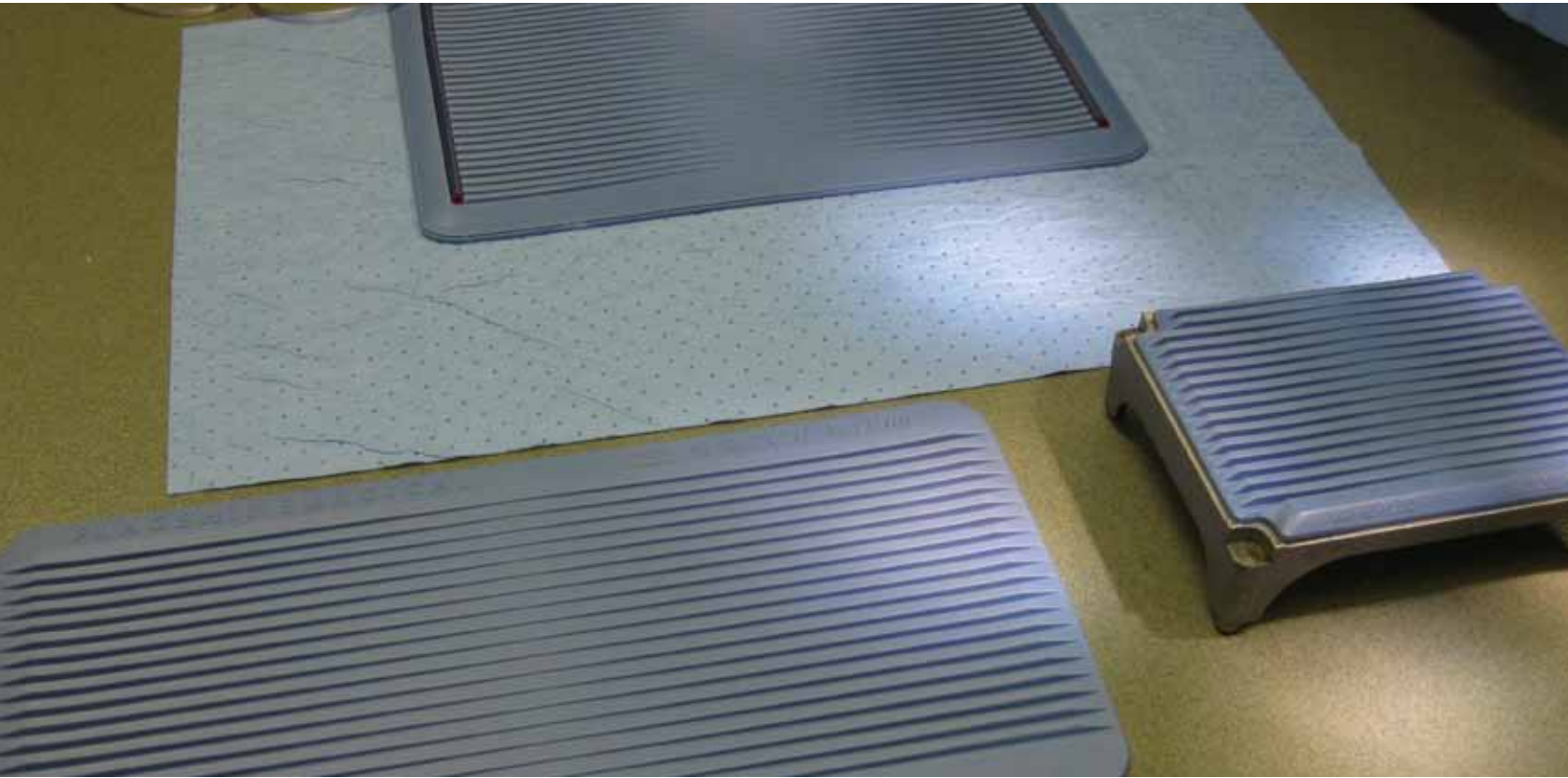

Suction Mats & Anti-Fatigue Mats

SPECIALTY FLUID CONTROL PRODUCTS

THE SUCTIONER™

- Anti-microbial disposable surgical mat
- Address three major OR issues: surgical comfort, fluid buildup, and OR Safety
- Designed and tested by world class surgeons, OR nurses, and technicians
- Incorporates revolutionary engineering process consisting of non-slip molded foam with custom tapered fluid channels
- Fluid channels help aid in keeping surgical team's feet dry and helps reduce the risk of tracking body fluids throughout the OR suite
- Elasticity and cushioning provide exceptional comfort during prolonged, strenuous surgeries; mitigating foot, leg, and back pain
- Incorporates a proprietary, highly-effective vacuum system for quiet, continuous suction. Attach directly to suction can or suction device using standard suction tubing
- Helps reduce need for countless blankets and towels placed on OR floor to absorb fluids
- Less mess on the OR floor helps provide shorter room turn-over time between surgeries
- Latex free



Patent Pending

Product #	Description	Quantity
USA-X10D	THE SUCTIONER™ (36" x 60")	10 mats/Cs

Ideal for the following surgeries:

- Arthroscopy
- Cardiothoracic Surgery
- Orthopaedic Surgery
- Trauma Surgery
- Transplant Surgery
- General Surgery
- Urologic Procedures
- Plastic Surgery
- Laparoscopic Surgery
- Cystoscopy
- Labor and Delivery & C-Sections

Your Front-Line Defense for Optimum Fluid Control

SPECIALTY FLUID CONTROL PRODUCTS



THE SURGICAL MAT™ and THE MINI MAT™

- Patented, anti-microbial disposable surgical mat
- Address three major OR issues: surgical comfort, fluid buildup, and OR Safety
- Designed and tested by world class surgeons, OR nurses, and technicians
- Incorporates revolutionary engineering process consisting of non-slip molded foam with custom tapered fluid channels
- Fluid channels help aid in keeping surgical team's feet dry and helps reduce the risk of tracking body fluids throughout the OR suite
- Elasticity and cushioning provide exceptional comfort during prolonged, strenuous surgeries; mitigating foot, leg, and back pain
- Latex free

Product #	Description	Quantity
USA-010D	THE SURGICAL MAT™ (20" x 39")	10 mats/Cs
USA-A10D	THE MINI MAT™ (13" x 17")	10 mats/Cs

Ideal for the following surgeries:


- Arthroscopy
- Laparoscopy
- Spine Surgery
- Total Joint Replacement & Revision Procedures
- Trauma Surgery
- Transplant Surgery
- Neurosurgery & Craniotomy
- Urologic Procedures
- Labor and Delivery & C-Sections
- Cardiothoracic Surgery
- Bariatric Surgery
- Plastic Surgery
- Oral Maxillofacial Procedures
- Ear, Nose, and Throat Surgery

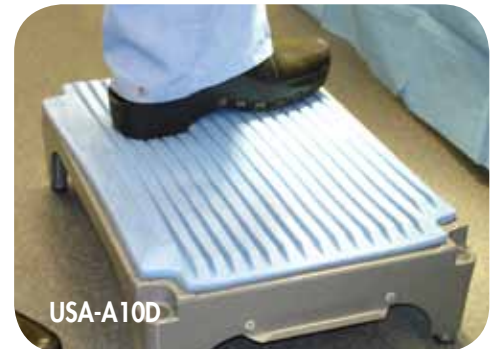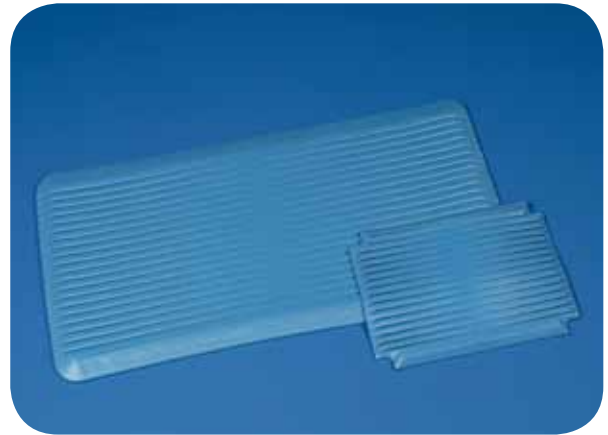

U.S. Patent No. 7,291,376

Suction Mats & Anti-Fatigue Mats

SPECIALTY FLUID CONTROL PRODUCTS

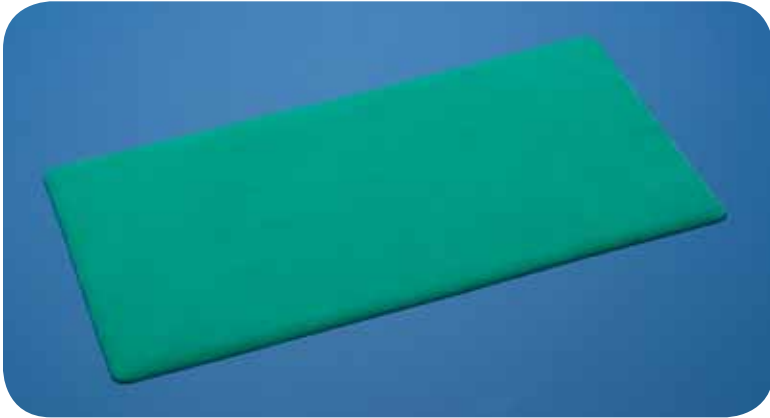

THE SMOOTH MAT™

- Disposable, economical, ergonomic Surgical Mat designed for use in the OR specifically for longer, dry surgical procedures
- Eliminates the exhaustive pain and discomfort from standing for extended periods of time in any OR
- Non-Sterile
- Latex free



Product #	Description	Quantity
USA-S10D	THE SMOOTH MAT™ (20" x 39")	10 mats/Cs

"The earliest and most common symptoms from prolonged standing during surgical procedures are discomfort and fatigue in the feet and legs - that is, in the shins, calves, knees and thighs. Pain and discomfort may also be felt in the hips, neck and lower back"...

"General recommendations: Anti-fatigue mats should be placed on the floors and standing stools"..."Facilities should also be aware of infection control issues for non-disposable anti-fatigue matting."

AORN Journal: Ergonomic Tool 4, June 2011

Your Front-Line Defense for Optimum Fluid Control

SPECIALTY FLUID CONTROL PRODUCTS



THE Rx1 MAT™

- Revolutionary anti-fatigue properties
- Helps reduce the strain, stress, and discomfort of prolonged standing in the workplace
- Both OSHA and AORN clearly advocates the use of anti-fatigue matting in the workplace to reduce musculoskeletal injury and improve workplace ergonomics
- Reusable - simply clean, as needed, with a mild detergent and water
- Non-Sterile
- Latex free

Product #	Description	Quantity
USA-R10D	THE Rx1 MAT™ (20" x 39")	10 mats/Cs

Ideal for the following work environments:

- Exam Rooms
- Surgical Facilities
- Work Stations
- Sterile Processing Departments
- Radiology Departments
- GI Suites
- Laboratories

Suction Mats & Anti-Fatigue Mats

SPECIALTY FLUID CONTROL PRODUCTS



Wheeled Storage Unit™

- Pre-Assembled unit with snap-on wheelcasters
- Convenient storage bin to stock surgical mats
- Holds 10 large mats and 8 mini mats
- Rollable cart to reduce lifting
- Latex free

Product #	Description	Quantity
USA-10WD	Wheeled Storage Unit™	1/Cs

Your Front-Line Defense for Optimum Fluid Control

CRITICAL 2-STEP DRAPING METHOD



Step 1: Following standard patient draping, twist the drapes around each other at the bottom, directly above THE SUCTIONER™.

- Fluids follow the pathway of the drapes. By twisting the drapes around each other, a V-shape or funnel-shape is created with the drapes above, which now aims downward toward THE SUCTIONER™. This will ensure that a majority of the fluids will hit the mat and be suctioned away.



Step 2: Take any standard clamp and secure the twisted drapes so that the clamp hangs over THE SUCTIONER™, ideally 6-12 inches above the mat.

- Note: Occasionally, it may be necessary to add additional drape clamps to redirect the fluids toward the mat during the surgical procedure, as the drapes may shift during surgery.

Suction Mats & Anti-Fatigue Mats



Improving Care. Improving Business.

DeRoyal Industries, Inc.
200 DeBusk Lane
Powell, TN 37849
888.938.7828 or 865.938.7828
www.deroyal.com

©2012 DeRoyal Industries, Inc.
Reprint # 0-2219
Rev. 1/12

The Suctioner, The Surgical Mat, The Mini Mat, Wheeled Storage Unit, The Smooth Mat, and The Rx1 are trademarks of Flagship Surgical, LLC.