


Lumitech Incisional Wound Dressing


The Lumitech multi-layer incisional wound dressing is intended to manage post-surgical wounds with low to medium exudate levels. The dressing incorporates **Lumina® Adhesive Technology**, a patented UV switchable adhesive that supports long wear time and a gentle dressing removal. The dressing is compatible with the Prospera Flex® NPWT System.

Innovative Features

- Peel and place, sterile multi-layer absorbent design
- Integrated Lumina Adhesive Technology
- Optimized design for application over various incision site locations
- NPWT compatible with dome site positioning flexibility
- Available in a variety of sizes

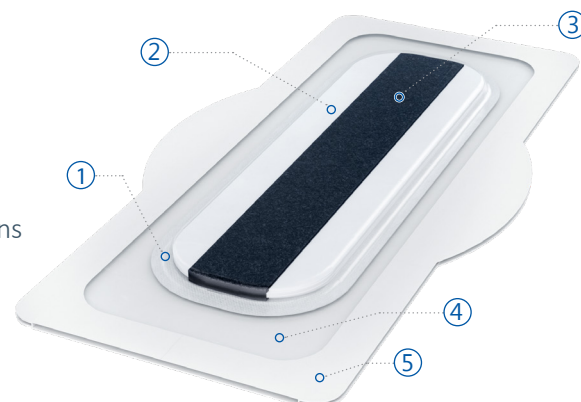

LUMITECH™

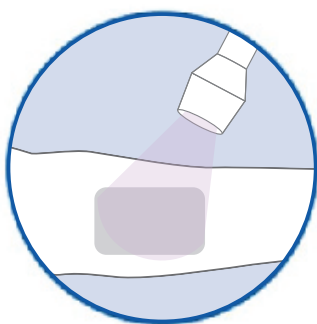


Multi-Layer Dressing Composition

1. Integrated wound contact layer with extended edges for increased skin integrity and patient comfort
2. Non-woven fibers support up to 20ccs of exudate absorption at the wound site when used independently or with NPWT¹
3. Polyurethane foam assists with pressure distribution when used with NPWT
4. Polyurethane film incorporating Lumina Adhesive Technology provides flexible, breathable, water-resistant, contaminant-blocking protection with a long wear time and gentle removal
5. Release liners and carrier layers for ease of application

Lumina Adhesive Technology Benefits²

- Adhesive adheres as strong as acrylic dressings, supporting strong seal in unswitched state
- Adhesive releases more gently than silicone following UV light exposure, minimizes pain and periwound trauma at dressing removal

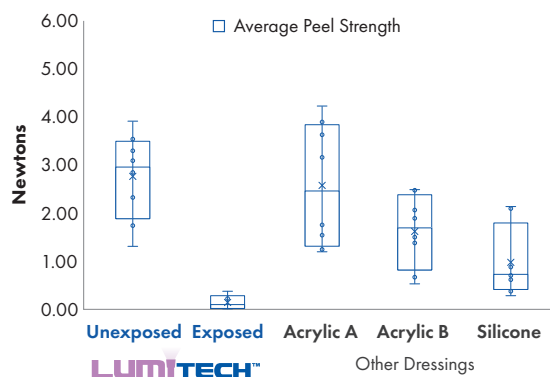

Compatible with the DeRoyal Prospera Flex NPWT Single-Use System



Product No.	Description	Packaging
46-LT20	Lumitech 20cm Incisional Dressing	10 Ea/Cs
46-LTUV	Lumitech Compatible UV Flashlight	10 Ea/Cs
46-LT20S-01-120	Lumitech Securement Strips (20cm)	6 Ea/Pouch, 20 Pks/Cs
46-LT20S-30	Lumitech Securement Strips (20cm)	6 Ea/Pouch, 5 Pks/Cs
47-820K-01	Incision Management Convenience Kit (Medium) x1 Prospera Flex® NPWT System, x2 Lumitech Incisional Wound Dressings x1 Lumitech Compatible UV Flashlight Accessories	1 Ea
47-820K	Incision Management Convenience Kit (Medium) x1 Prospera Flex NPWT System, x2 Lumitech Incisional Wound Dressings x1 Lumitech Compatible UV Flashlight Accessories	5 Ea/Cs

**To learn more about our products or to place an order,
Call 1-888-938-7828 or visit [DeRoyal.com](https://www.DeRoyal.com).**

1- Data on file, 20cm Dressing

2- Cayce, Jonathan, Ph.D., MS, Goins, Lily, & Warrington, Gavin. (2023, April 26). Light-Deactivated Adhesive: The Solution to Maintaining Dressing Integrity and Protecting Skin. Society for Advanced Wound Care, National Harbor, MD.

DeRoyal, the DeRoyal logo, Prospera Flex, and Lumitech are trademarks or registered trademarks of DeRoyal Industries, Inc. Lumina is a registered trademark of Lumina Adhesives AB.